

bebird pro

⚠ Note: Do not use the charging equipment with output voltage exceeding 5V

Step 1 : Download the APP

iOS : Search "bebird" in App Store

Android: Search "bebird" in Google Play




Or scan the QR code to download the app



Step 2 : Connect to the WiFi

1 Open APP Connect to the WiFi in the settings

2 Press the power button to turn on

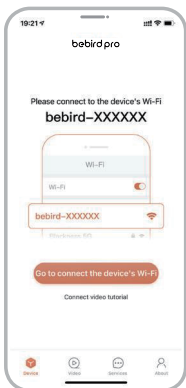

— The blue light
flashes after
starting


3

Connect Product Wi-Fi
"bebird-R1-XXXXXX"



4

Back to bebird app
start to use it

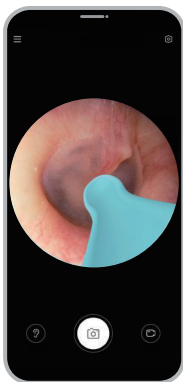

Accessory replacement

Please check that
the ear tip is tight
before using it.



The ear tip has its
own thread and
can be unscrewed
clockwise for
replacement.

For a better user experience, we recommend
replacing the ear tip at least every 3 months .
(This can be adjusted according to the frequency of use.)

Accessory

1

PC ear spoon
x 1



2

Cotton swab
ear tip x 1



3

Charging
Cable x 1



4

Storage tube
x 1



Troubleshooting

Problem	Solution
No response after startup	Please use the incidental charging cable to charge. The charging time is 50mins
The app exits quickly	Delete the existing app and scan QR to reinstall it
Feel tricky with replacing ear tips	Align the replacement ear tip to camera and tighten clockwise
Indicator light flashes	Please use the incidental charging cable to charge. The charging time is 50 mins
Abnormal display in app	Interrupt the connection, restart the product to connect again
If you still have problems with the above procedures, please contact our Customer Service Representatives directly for further assistance	

Hazardous substance

Part name	Hazardous material					
	(Pb)	(Hg)	(Cd)	(Cr(VI))	(PBB)	(PBDE)
Printed Circuit Components	×	○	○	○	○	○
Case	○	○	○	○	○	○
Plastic rubber accessory	○	○	○	○	○	○

This form is formulated in accordance with SJ/T 11364

○ : represent the hazardous material ratio is below GB/T26572-2011

× : represent the hazardous material ratio is over GB/T26572-2011

Product parameters

Product Name	Ear Wax Removal Tool
Pixel	8 Mega-pixels
Depth of field	10-50mm
Lens diameter	3.5mm
Image transmission rate	20fps
Battery	180mAh lithium battery
Battery life	70 minutes
Charging time	50 minutes
Input current	DC 5V 0.5A
Gyroscope	6 axis
Network standard	IEEE 802.11b/g
Working frequency	2.4GHz
Rod weight	13g
Image sensor	CMOS
Working temperature	-10 ~ 50°C
Optimum focal length	13.5mm



Attention

- 1、 When cleaning lens, please wipe it carefully with professional alcohol cotton swabs.
- 2、 Please observe your surroundings before use. Stay at least 2ft away from people, animals, or moving objects when using . Don't use in a moving vehicle or other unstable environments. Please don't move around when using.
- 3、 This product is unsuitable for children under 3 years old.
- 4、 This product is not suitable for children to use alone.
- 5、 Do not place in a humid environment for a long time or soak in liquids.
- 6、 There is a built-in rechargeable lithium-ion battery. Even if you do not use it, you should charge it once a month to avoid long-term battery damage.
- 7、 Please use a 5V charger to charge the device. Do not use this product while being charged.
- 8、 Avoid direct exposure to the sun.
- 9、 Do not put heavy objects on the product.
- 10、 It is normal for the temperature to rise slightly when using (up to 95°F).

Warranty

User info	
User name:	Telephone number:
Address:	
Seller:	Telephone number: