

bebird pro

Ear-care
You Can See



Please consult the instructions before using this product

Contents

1. English version instruction manual	EN-1 --- EN-5
App Download	EN-1
App Connection	EN-1
Device Overview	EN-2
Accessory Installation	EN-3
Bebird EarSight Plus Series Accessory List	EN-4
Troubleshooting	EN-5
Product Specifications	EN-5
WARNINGS	EN-5

App Download



or



Search for "**bebird**"
in your phone's app store.

Scan the QR code to download
and install the **bebird** app.


App Connection

- 1 Press and hold the power button to turn on the device.




Once the device powers on,
a blue light will illuminate.


- 2 Open the Bebird App, and click "Connect the device"

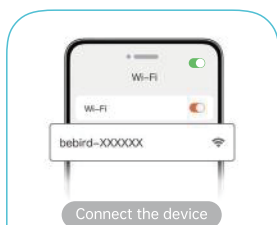

- 3 Connect Product Wi-Fi
"bebird-xx-xxxxxx"

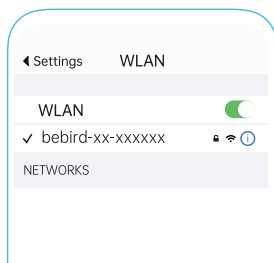

- 4 Return to the Bebird app and click "Begin to Clean Ears" to start using the device.

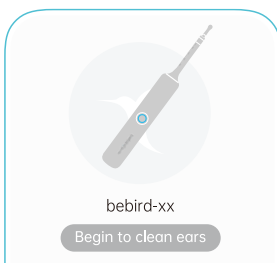

Device Overview

Silicone
Covered
Tips

Camera Lens


Hose

Attention:

Please note that our hose is designed with a specific bending direction. When using, ensure the scale side is facing forward.

The hose should only be bent **left** and **right** to avoid damage or reduced functionality.

Do not bend in other directions, otherwise it may cause damage to the product.

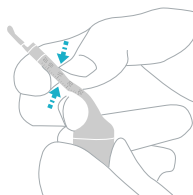

Do not use brute force. Before using the tweezers and ENT set, the bent lens needs to be straightened again.

Power Button

When the device is powered on, a blue light will illuminate.


☀ During power-on and connection, the blue indicator light will flash.

● Once connected, the blue light will remain lit.

● While charging, the indicator light will glow red.

○ When charging is complete, the red indicator will turn off.

☀ The blue indicator will flash rapidly when the battery is low.

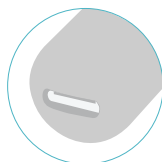

Charging Port

Accessory Installation

Protective Plug Installation:

Place the protective ear canal plug securely over the camera hose, ensuring it is positioned correctly.

This will prevent the device from penetrating too deeply into the ear canal.



Ear Spoon Installation:

To install the ear spoon, simply screw the metal-covered ear spoon tightly in place. This method applies to all types of ear spoons, including transparent two-color ear spoons, transparent ear spoons, ring-shaped rubberized ear spoons, and sponge ear spoons.




Tweezer Mounting:


Place the tweezers extension onto the body of the device. Tightly screw the head into place to secure.



Bebird EarSight Plus Series Accessory List

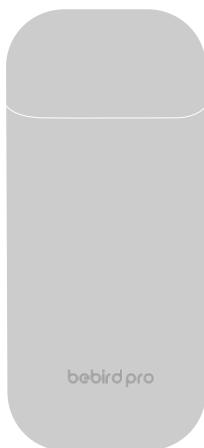

Accessory Case

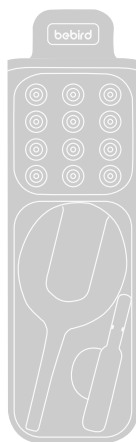

Accessory Tray



Nasal Speculum



Ear Canal Speculum



Tongue Depressor



Metal, Silicone Wrapped
Ear Spoon



Double-color Q-elastic
Ear Spoon



Double-color Annular
Ear Spoon




Sponge Ear Spoon



Plastic Ear Spoon



Tweezers
Component-B

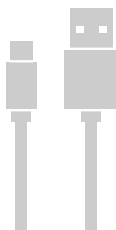

Charging Cable



Ear canal
Protective Plug



Lens Brush



Tweezers
Component-A

Troubleshooting

Problem	Solution
No response after startup	• Use the enclosed charging cable to charge
The LED light of the camera is not on	• Interrupt the connection, restart the device to connect again
The app closes quickly	• Delete the existing app and reinstall
The indicator light blinks quickly	• Use the enclosed charging cable to charge
Abnormal display in app	• Interrupt the connection, restart the device to connect again
The wifi connection was unstable and disconnected after a few seconds	• Interrupt the connection, restart the device to connect again
Image is blurry	• Wipe the camera lens with a double-ended cleaning brush, cleaning cloth, or alcohol swab
Movement feels unnatural	• Switch the left/right ear mode in the app interface
If the issue persists after following the above operations, please contact our customer service, and we will provide you with further effective solutions.	

Product Specifications

Brand: Bebird Pro

Name: Ear Wax Removal Tool

Usage time: Up to 80 minutes

Charge time: About 50 minutes

Intelligent orientation: Gimbaled Intelligent Gyroscope

Image sensor: CMOS

Input Current: 5V $\equiv$ 0.5A

Operating Frequency: 2.4GHz

Lens grade: 3P* Slim Lens, IP67 Waterproof, Dustproof & Anti-Fog

Operating Environment: -10°C to 50°C

Image Frame Rate: 20fps

Network Standard: IEEE 802.11b/g

Optimum Focal Length: 12.5mm



WARNINGS

- **Cleaning:** Wipe the lens with a lens brush.
- **Assembly:** Ensure all device accessories are securely attached to prevent them from falling off during use.
- **Operating Environment:** Use the device in a suitable, private environment. Avoid using it in crowded spaces.
- **Battery Maintenance:** The device includes a built-in rechargeable battery. If unused for extended periods, charge it once a month to maintain battery health.
- **Supervision for Children:** Children should not use this device unsupervised to avoid accidental injury.
- **Sunlight Exposure:** Keep the main device and accessories out of direct sunlight to prevent softening or damage.
- **Water Exposure:** Do not submerge the device in any liquids.
- **Temperature Notice:** During use, the device may warm slightly due to the LED lights; this is normal and is capped at approximately 45°C (113°F).
- **Intended Use for Tweezers:** Use the tweezers only to remove earwax or other foreign objects from the ear canal.

