

**bebird pro** *uena*

See the best in you



Please refer to the instructions before using this product

# Contents

|                                             |               |
|---------------------------------------------|---------------|
| 1. English version instruction manual ..... | EN-1 --- EN-9 |
|---------------------------------------------|---------------|

## Download App



or



**iOS** Search "**Bebird**"  
in App store

Scan QR code to  
download the App

Watch Tutorial Video  
via this QR code

**Android** Search "**Bebird**"  
in Google Play store

## Auto Connection

1


Turn on the device and  
open the mobile App.

2

Click "Connect device"  
on the homepage.

3

After the search is  
successful, click "Connect".

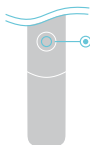

Press and hold  
for 2s to power  
on; the blue  
light will flash.

### < Connect device

Auto search

Connect manually

### < Connect device

1 bebird-xx-xxxxxx 2



Connect

## Manual Connection

1 Click "Connect Manually" on  
the "Connect Device" page.

2 Go to your phone's Settings > Wi-Fi,  
and connect your device's Wi-Fi.

### < Connect device

Auto search

Connect manually

### < Wi-Fi

Wi-Fi



✓ bebird-xx-xxxxxx



NETWORKS

## Device Overview

Silicone Covered Tips

Camera Lens

Hose

bebird pro

### Power Button

When the device is powered on, a blue light will illuminate.

- ☀ During power-on and connection, the blue indicator light will flash.
- Once connected, the blue light will remain lit.
- While charging, the indicator light will glow red.
- When charging is complete, the red indicator will turn off.
- 🔋 The blue indicator will flash rapidly when the battery is low.

Charging Port

## Accessory Installation

### Protective Plug :

Place the ear canal protective plug in the appropriate position on the hose to prevent it from being inserted too deeply into the ear canal.



### Ear Spoon:


Unscrew the current accessory by twisting it, counterclockwise to loosen and remove it.

Attach the desired accessory by screwing it in clockwise until its tight and secure.




## How to Use the Tweezers Components?

Place tweezers component-A onto the device body firstly, then install and tighten tweezers component-B

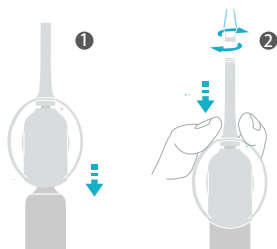

Gently pinch the handle of the tweezers with two fingers to close the tweezers.

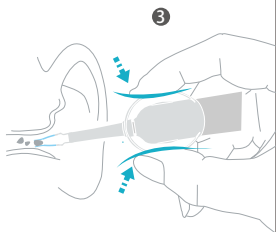

### Warm Reminder

Install the tweezers component correctly before use to avoid its loosening or falling off.


## How to Use the Hose lens?

In most cases, you don't need to bend the hose to use it. If necessary, please note that:

Facing the scale, slowly bend the hose with your fingers to the left and right sides of the scale.


Do not bend in other directions, otherwise it may cause damage to the product.

Do not use brute force.  
Before using the tweezers, the bent hose needs to be straightened again.

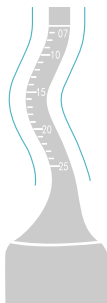

## How to Use the Tacky Head or Sponge Head?

Attach the Tacky Head or Sponge Head to the hose securely, prevent loosening. Ear debris can be removed by Stubbord Ear debris can be removed by Tacky Head. Oily earwax can be absorbed by the Sponge Head.

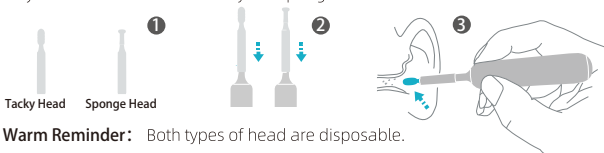


**Warm Reminder:** Both types of head are disposable.

## How to Use the Dropper?

- 1 Insert the nozzle into the liquid and gently press the air bag to suck up the liquid.
- 2 Install the dropper on the device securely




- 3 With the help of a camera, liquid can be precisely dripped into the ear.

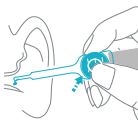

### Warm Reminder :

1. Please avoid immersing the lens in liquid, otherwise it may cause the lens to become blurred and need to be wiped before continuing to use.
2. This dropper is for use with ear drops only (this product does not contain ear drops). If used with other liquids, please seek professional advice.

## Accessory List 1

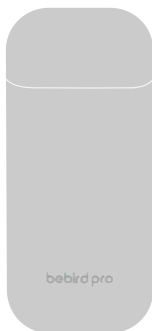

**Travel Box**



**Main Device**



**Inner Tray**



**Tacky Head**



**Smooth round-tipped metal ear spoon**



**Dome Ear Spoon  
(for practice use)**



**Plastic Ear Tip  
(hard)**



**Metal, Silicone  
Wrapped Ear Tip  
(Moderate)**




**Double-color  
Q-elastic Ear Tip  
(soft)**



**Sponge Head**

## Accessory List 2

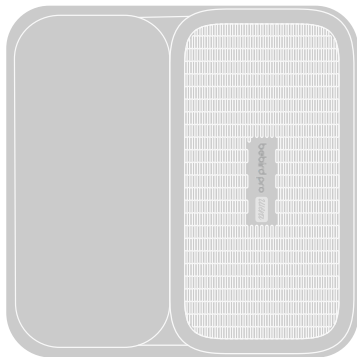

**Master Box**



**Type-C  
Charging Cable**



**Ear Canal  
Speculum**



**Protective  
Plug**



**Lens Brush**



**Tweezers  
Component-B**



**Nasal  
Speculum**



**Tongue  
Depressor**




**Dropper**



**Tweezers  
Component-A**

## How to Use the Ear/Nasal/Throat Observation tools?


Put on the ear canal speculum tightly, and slowly go into the ear.

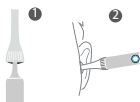

Put on the tongue depressor tightly to prevent it from loosening.



Put on the nasal speculum tightly, and slowly go into the nose.

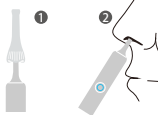

**\*Note:** Please make sure the tool is securely attached. If feel uncomfortable, please do not continue.

## Troubleshooting

| Symptoms                                                                         | Solution                                                                                                                                                                                                                                 |
|----------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Camera light off                                                                 | Check if the brightness setting in the settings is set to the lowest setting.                                                                                                                                                            |
| Lens blurry                                                                      | Gently wipe the lens surface with the included lense brush.                                                                                                                                                                              |
| App crashes                                                                      | Uninstall the existing app and reinstall.                                                                                                                                                                                                |
| No response on startup                                                           | Please charge it with the included charging cable.                                                                                                                                                                                       |
| Indicator light flashes rapidly                                                  | Battery is low, please recharge the battery in time.                                                                                                                                                                                     |
| The product fails to display live video after successfully connecting to the app | Check if the app permissions are disabled.                                                                                                                                                                                               |
| Wi-Fi connection is unstable, sometimes disconnecting within seconds             | Restart the device and reconnect to the app.                                                                                                                                                                                             |
| Indicator light and lens light flash simultaneously                              | Flashing indicates a firmware upgrade to be done. Firmware upgrade will take about 30mins, please charge the device for 30 mins to avoid power cut during the upgrade. Once upgrade is done, no more flashing, and ready for normal use. |
| Incorrect orientation during use                                                 | Switch left/right ear mode on the bottom of the app.                                                                                                                                                                                     |

If the issue persists after following the above operations, please contact our customer service, and we will provide you with further effective solutions.

## Product Specifications

Brand: **Bebird Pro**

Model: **EarSight Ultra**

Intelligent orientation: **Gimbaled**

Intelligent Gyroscope

Usage time: **Up to 80 minutes**

Network Standard: **IEEE 802.11b/g**

Image sensor: **CMOS**

Operating Environment: **-10°C to 50°C**

Input Current: **5V  $\approx$  0.5A**

Operating Frequency: **2.4GHz/5GHz**

Optimum Focal Length: **12.5mm**

Image Frame Rate:  **$\geq$  25fps**

Lens grade: **4P\* Slim Lens, IP67 Waterproof**



## Warnings



- **Cleaning:** Wipe the lens with a lens brush.
- **Assembly:** Ensure all device accessories are securely attached to prevent them from falling off during use.
- **Operating Environment:** Use the device in a suitable, private environment. Avoid using it in crowded spaces.
- **Battery Maintenance:** The device includes a built-in rechargeable battery. If unused for extended periods, charge it once a month to maintain battery health.
- **Supervision for Children:** Children should not use this device unsupervised to avoid accidental injury.
- **Sunlight Exposure:** Keep the main device and accessories out of direct sunlight to prevent softening or damage.
- **Water Exposure:** Do not submerge the device in any liquids.
- **Temperature Notice:** During use, the device may warm slightly due to the LED lights; this is normal and is capped at approximately 45°C (113°F).
- **Intended Use for Tweezers:** Use the tweezers only to remove earwax or other foreign objects from the ear canal.
- The inside of the ear canal is relatively fragile. To avoid damaging the ear canal, please do not use this product to scrape the ear canal vigorously.

## Warranty Service

### Customer information

Customer Name:

Customer Telephone:

---

Retailer Name:

Customer Address:

---